

Corning® Supersomes™ Ultra HEK293-derived Human Aldehyde Oxidase (AOX1) Cytosol

Catalog Number.....456801
 Lot Number.....6342021

Storage Conditions.. -80°C
 Date Released.....2016 December
 Expiration Date.....2021 December

Product Description: Human AOX1 cytosol is prepared from Human Embryonic Kidney (HEK) 293 cells overexpressing human Aldehyde Oxidase (AOX1) (Gene Accession No. NM_001159). Aldehyde oxidase is a cytosolic enzyme that plays an important role in xenobiotic biotransformation. Primarily expressed in liver, AOX1 has broad substrate specificity, catalyzing the oxidation of aldehydes, aromatic azaheterocycles and nitrogen-carbon double bond containing compounds.

Negative Control Cytosol (Cat. No. 456800) prepared from HEK293 cells transfected with empty vector can be used as a control for this product.

Package Contents.....2.5 mg protein in 0.5 mL
 Protein Content.....5.0 mg/mL in 0.1M potassium phosphate buffer, pH7.4

QUALITY CONTROL:

Specification	Criteria	Result
Enzymatic Assay	Ultra Supersomes: Human AOX1	PASS
	Substrate: 250 µM zaleplon	
	Incubation time: 60 min	
	Activity (pmol/mg/min) >=2 pmol/mg/min.	

METHOD (zaleplon assay): The reaction is initiated by adding 10 µL of 5 mg/mL AOX1 cytosol into 90 µL 100 mM potassium phosphate (pH 7.4, 37°C) containing 278 µM zaleplon substrate. The reaction mixture is incubated for 60 min at 37°C. The reaction is terminated by adding 100 µL pre-chilled acetonitrile containing 0.4 µM labetalol as internal standard. The metabolite formation is quantified using LC-MS analysis.

STORAGE and HANDLING

- Thaw on ice and keep on ice until use.
- The activity can be compromised upon freezing and thawing. The activity remains 80% or above if the sample is thawed less than 4 times. After initial thaw, it is highly recommended to aliquot sample into single use vials to minimize freeze-thaw cycles.
- Triplicate samples per assay conditions are recommended.
- The concentration of test compound and incubation time usually depends upon the specific characteristics of the drug.

For Research Use Only. Not for use in diagnostic or therapeutic procedures.

For a listing of trademarks, visit www.corning.com/lifesciences/trademarks

© 2017 Corning Incorporated

Discovery Labware, Inc.
6 Henshaw Street
Woburn, MA 01801
Tel: 1.978.442.2200 (U.S.)
CLSTechServ@Corning.com

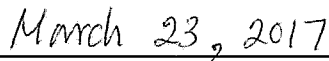
CORNING

www.corning.com/lifesciences

Quality Assurance:

A handwritten signature in black ink, appearing to read "Hundo", is written over a horizontal line.

Date:

The date "March 23, 2017" is handwritten in black ink over a horizontal line.

SAFETY RECOMMENDATIONS:

Safety assessment indicates this product is not hazardous, therefore no SDS (Safety Data Sheet) is provided. Handle in accordance with good industrial hygiene and laboratory safety practices.

For Research Use Only. Not for use in diagnostic or therapeutic procedures.

For a listing of trademarks, visit www.corning.com/lifesciences/trademarks
© 2017 Corning Incorporated